

LiNA SeaStar™ and LiNA Skin Hook™

Complete single-use system for
retracting soft tissue from the
operative surface




Elastic Retraction System



Innovation in Gynecology

LiNA SeaStar™ and LiNA Skin Hook™

The LiNA SeaStar™ and LiNA Skin Hook™ offer a single-use solution for retracting delicate soft tissue. The flexible design of the system allows adjustment and repositioning of the LiNA SeaStar™ during the surgical procedures and provides exposure of the surgical site.

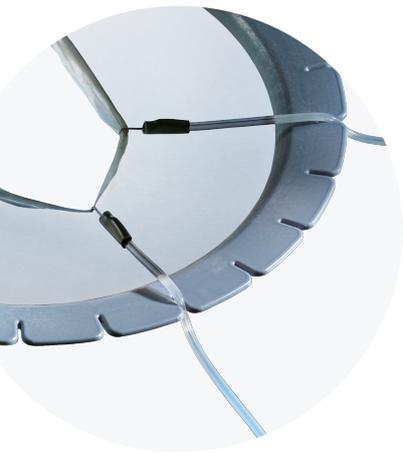

Visibility. Self-retaining.

Freeing the hands of the surgeon

The LiNA SeaStar™ and LiNA Skin Hook™ provide improved surgical space and visibility. Allowing you to focus on operating, teaching and training rather than retraction.


Versatility. Broad range of available options.

LiNA SeaStar™


The LiNA SeaStar™ is provided in a broad range of shapes and sizes. The frames can be adjusted to fit various anatomical areas, thereby facilitating different surgical applications e.g. vaginal and rectal surgery and other surgical disciplines.

LiNA Skin Hook™


The LiNA Skin Hook™ is available with three different stainless steel hook options: small sharp, small blunt, large blunt. The flexible, elastic bands can be adjusted during the procedure to create optimal tension.

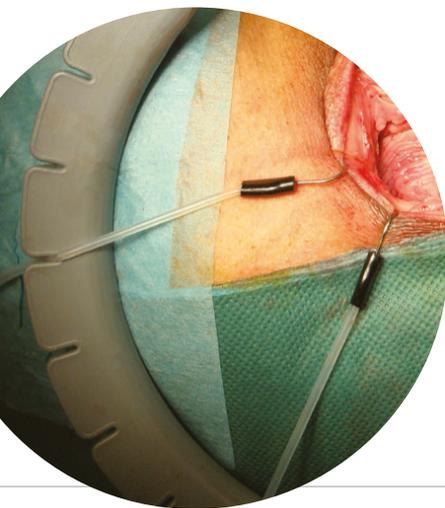


Prepare and drape the surgical area in usual manner. Apply the LiNA SeaStar™ according to the following steps:

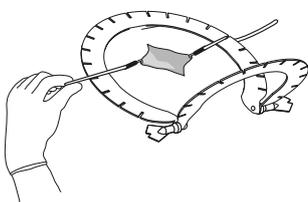

1. Place the LiNA SeaStar™ around the initial incision and insert the LiNA Skin Hook™ in directly opposing quadrants with equal tension.

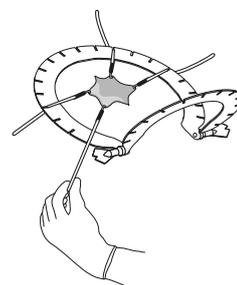

2. As dissection progresses, move the position of the LiNA Skin Hook™ into the deeper layers of fascia. You may then add new LiNA Skin Hooks or reposition those already in place to achieve full balanced retraction from any angle.

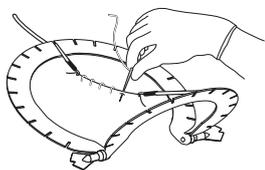


3. To ease closing, maintain constant opposing tension.




LiNA SeaStar™ and LiNA Skin Hook™

The LiNA SeaStar™ and LiNA Skin Hook™ are single-use.

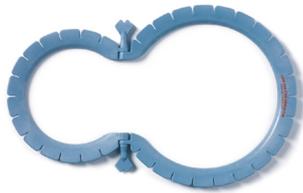

REF 4000-S-SL

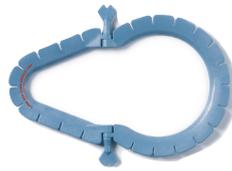


REF 4000-S-VS



REF 4000-S-SS



REF 4000-S-VV

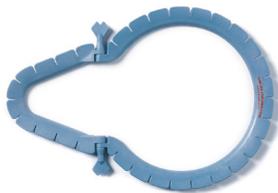

REF 4000-S-VL



REF 4000-S-Q

LiNA Skin Hooks
are actual size



REF SH-110
Small sharp hook



REF SH-120
Small blunt hook



REF SH-220
Large blunt hook

Ordering information:

LiNA SeaStar™ and LiNA Skin Hook™



Product Code	Description	UOM
4000-S-SL	LiNA SeaStar™ – Small/large. Single use.	Box of 5, sterile
4000-S-SS	LiNA SeaStar™ – Small/small. Single use.	Box of 5, sterile
4000-S-VL	LiNA SeaStar™ – V-shape/large. Single use.	Box of 5, sterile
4000-S-VS	LiNA SeaStar™ – V-shape/small. Single use.	Box of 5, sterile
4000-S-VV	LiNA SeaStar™ – V-shape/v-shape. Single use.	Box of 5, sterile
4000-S-Q	LiNA SeaStar™ – Squared. Single use.	Box of 5, sterile
SH-110	LiNA Skin Hook™ – Small sharp hook. Single use.	2 pieces in a pouch, box of 5 pouches, sterile
SH-120	LiNA Skin Hook™ – Small blunt hook. Single use.	2 pieces in a pouch, box of 5 pouches, sterile
SH-220	LiNA Skin Hook™ – Large blunt hook. Single use.	2 pieces in a pouch, box of 5 pouches, sterile

Distributor:

©2019 LiNA Medical ApS. All rights reserved.

LiNA is a registered trademark of LiNA Medical in EU, US and CN.

LiNA Medical ApS

• Formervangen 5
DK-2600 Glostrup
Denmark

• Tel.: +45 43 29 66 66
Fax: +45 43 29 66 99
info@lina-medical.com

• www.lina-medical.com
CVR: 10054974

EN ISO 13485:2016
Quality Management
System

